



Syringes

SOL-CARE™ Safety Syringe



SOL-CARE™ Safety Syringe is a retracting needle manual safety syringe. It is designed to prevent needlestick injuries and syringe reuse by retracting the needle into the barrel.

- Compatible with all standard hypodermic needles.
- Easy to use.
- Fingers remain behind the needle during activation of safety mechanism.
- Clear and bold scale markings.
- Visual and audible confirmation of safety mechanism activation.

CE-marked



Product Options	Available Syringe Sizes
with Fixed Needle	0.5 ml 1 ml
with Exchangeable Needle	3 ml 5 ml 10 ml 20 ml
without Needle	3 ml 5 ml 10 ml 20 ml

Flexibility in product customization for packaging and labeling.

SOL-CARE™ Safety Syringe



Syringes

Quick and easy activation of safety mechanism.



PUSH



PULL



SNAP

Product Specification

CE-marked

SOL-CARE™ Safety Syringe product line offers comfortable and convenient injection solutions for key therapeutical areas.

Therapeutical area	Product	Syringe size
Diabetes	SOL-CARE™ Insulin Safety Syringe with Fixed Needle	0.5 ml 1 ml
Allergy	SOL-CARE™ TB Safety Syringe with Fixed Needle	1 ml
Immunization	SOL-CARE™ Safety Syringe with 26G 3/8" and intradermal bevel	1 ml

Flexibility in product customization for packaging and labeling.



Syringes

SOL-CARE™ Safety Syringe


SOL-CARE™ TB Safety Syringe with Fixed Needle

Reference number	Capacity	Gauge size	Outer diameter (mm)	Length		Shelf Box / Shipping Case (pcs)
				(mm)	(inch)	
100067IM	1 ml	28G	0.35 mm	12.5 mm	1/2"	100/1000
100019IM	1 ml	27G	0.4 mm	12.5 mm	1/2"	100/1000
100070IM	1 ml	26G	0.45 mm	10 mm	3/8" IDB	100/1000
100018IM	1 ml	25G	0.5 mm	16 mm	5/8"	100/1000



SOL-CARE™ Luer Lock Safety Syringe without Needle

Reference number	Capacity	Shelf Box / Shipping Case (pcs)
120006IM	3 ml	100/800
120007IM	5 ml	100/800
120008IM	10 ml	100/800
120009IM	20 ml	50/600

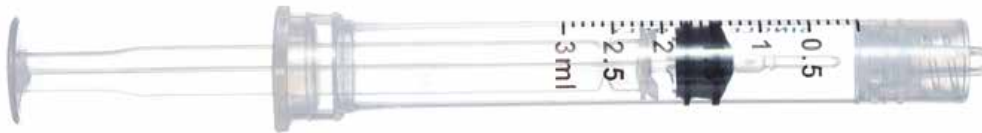

SOL-CARE™ Insulin Safety Syringe U-100 Insulin Only (with Fixed Needle)

Reference number	Capacity	Gauge size	Outer diameter (mm)	Length		Shelf Box / Shipping Case (pcs)
				(mm)	(inch)	
100089IM	0.5 ml	30G	0.3 mm	8 mm	5/16"	100/1000
100002IM	0.5 ml	29G	0.33 mm	12.5 mm	1/2"	100/1000
100017IM	1 ml	29G	0.33 mm	12.5 mm	1/2"	100/1000



More sizes available on demand